www.ocublink.com
chau@ocublink.com

SETUP GUIDE

OcuBall FBR

Thank you for choosing OcuBall for foreign body removal training. This guide shows you how to mount the support stand on your slit lamp and attach an OcuBall. Questions or suggestions? Email us any time.

1 Mount the support stand

Our updated stand is designed to fit most slit lamp head rests on the market.



Head of the stand
Pin that holds the OcuBall

Foot of the stand
Holes that line up with the pins on your chin rest

Figure 1 The OcuBall support stand

1. Place the stand on the head rest of your slit lamp (Figure 1).
2. Line up the holes in the foot of the stand with the holes in the chin rest.
3. Push the pins of your chin rest through each hole to secure the stand.
4. You are all set (Figure 2).



Figure 2 The support stand mounted on two different head rests

2 Attach the OcuBall

The OcuBall attaches to the pin at the head of the stand.

1. Take an OcuBall out of its sealed plastic container.
2. Line up the hole on the back of the OcuBall with the pin at the head of the stand, and push it on until snug.
3. Bring the OcuBall into focus and adjust the lighting of your slit lamp.
4. You are ready to go.

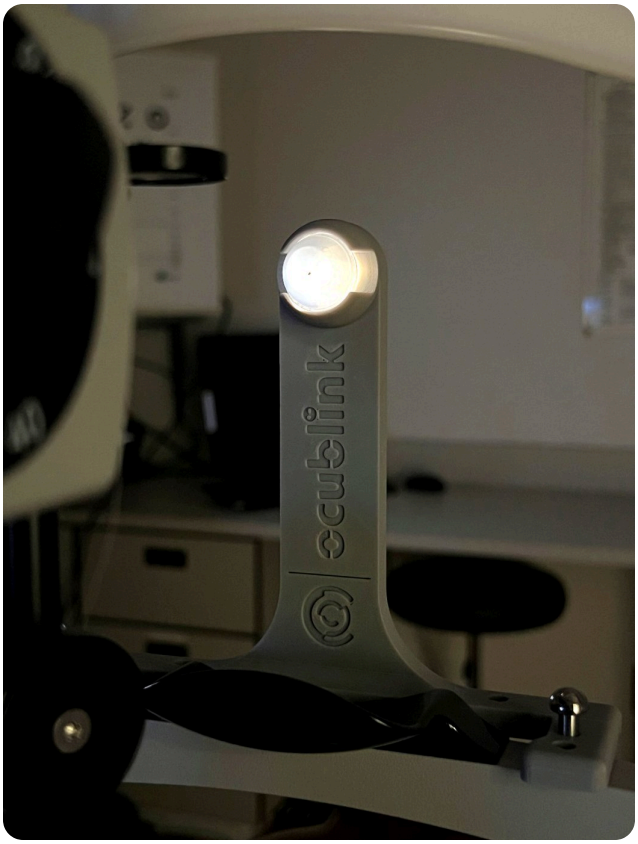


Figure 3 OcuBall set up on the slit lamp

Storage. Keep OcuBalls in the freezer for up to 180 days. Take out as many as you need 24 hours before use, so the eyes soften and the foreign bodies have time to rust. Each OcuBall is single use; the stand and face mask are reusable.

Happy foreign body removal!

Your OcuBlink team